



Xtraordinary - packaged as a sterile kit in three time-saving configurations, ready for all long bone trauma applications.

Xtra light - as strong as conventional fixators, only lighter, for easier application and more patient comfort.

X-ray friendly - the first large fixator body to be radiolucent.

XCaliber - High Caliber Fixation from Orthofix.



XCaliber. **X**actly.

- Surprisingly Light Weight
- Excellent Stability
- Sterile Kit, Ready to Use



XCaliber Hybrid


XCaliber Hybrid

Excellent periarticular fixation in small fragments and osteoporotic bone

Up to four wires on one ring for increased stability.



Rings convert from 2/3 to full circumference.

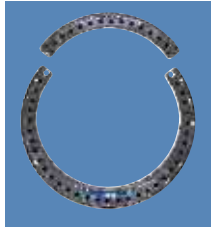

Radiolucent ring, fixator and stabilizing bars.



Ball-joints allow correction up to 44° from straight position, in any direction, with fixator in place.



Telescoping body for fracture reduction, and compression/distraction.



Ball-joints lock/unlock from either side. One instrument for fixator locking/unlocking.

Case Study

25-year-old male, metaphyseal fracture (AO type C2).



Minimal internal fixation of the articular surface with two cannulated screws and application of the XCaliber Hybrid fixator, with four tensioned wires proximally and two OsteoTite screws distally.



Fixator removed at 110 days.



Case provided by: Dr. M. Manca, Hospital of Massa-Carrara, Italy.



Pre-assembled and Packaged Sterile
Ready to use, straight from the box.

XCaliber Fixator Is Available In Three Configurations.

Use with XCaliber Bone Screws.



Meta-Diaphyseal



Hybrid



Articulated Ankle

Ordering Information

Part No.	Product Description
99-91080	XCaliber Hybrid Kit, Sterile
99-91038	Supplementary Screw Holder Kit, Sterile

Your Distributor is:

www.orthofix.com

XC-0309(A)-PL-E0 05/04

ORTHOFIX
ALWAYS INNOVATING