
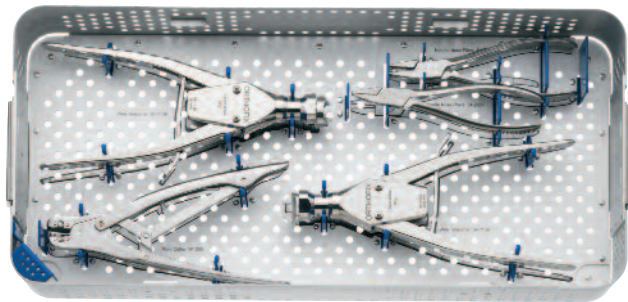


TRAYS


Instruments

Part Number
450178 C

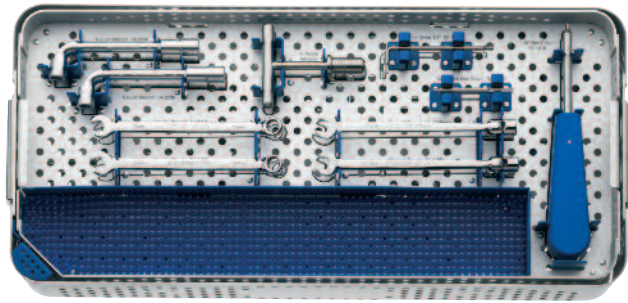


External Supports

Part Number
450179 C

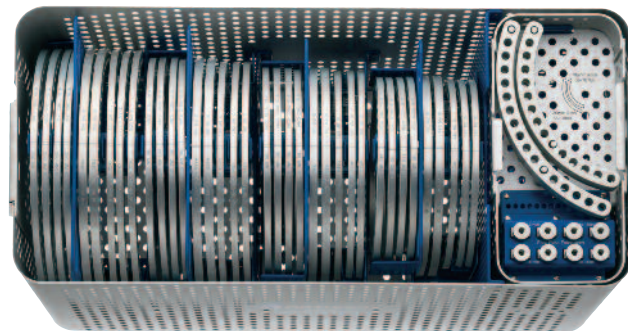


Fixation Elements

Part Number
450180 C

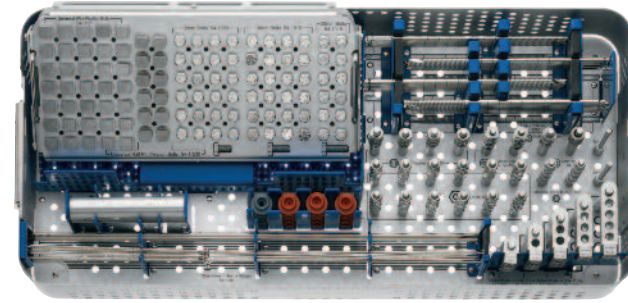

Hinges and Distractors

Part Number
450181 C

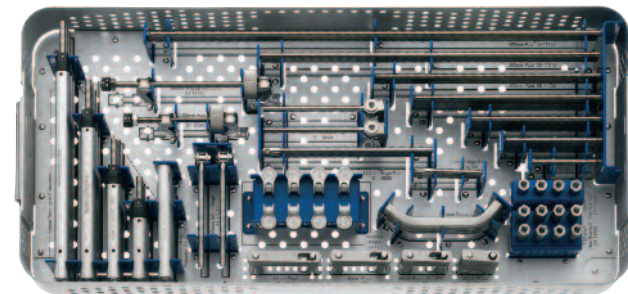

ORDERING INFORMATION

For specific part numbers and descriptions, refer to the TrueLok Parts Reference Guide, available upon request.



1.800.266.3349 (BoneFix)
www.orthofix.com

Adult & Pediatric Deformity Correction | Foot & Ankle | Bone Growth Stimulation | Fixation | Biologics


ORTHOFIX®

TL-0919 A-PL-US © Orthofix Inc. 11/2009

 DEFORMITY CORRECTION


TRUELOK™
RING FIXATION SYSTEM

TrueLok™ is a surgeon designed, user-friendly ring fixation system for limb lengthening, deformity correction and foot and ankle applications.



TrueLok is Simple
TrueLok is Stable
TrueLok is Versatile


ORTHOFIX®




TrueLok IS SIMPLE



Angular Distractor



Inline Hinge

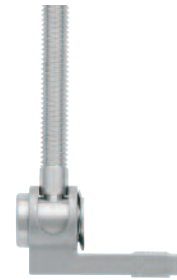

Outboard Hinge

TrueLok offers pre-assembled functional components which are easy to connect and operate:

- TrueLok inline and outboard hinges are dynamic connectors for deformity correction.
- A TrueLok angular distractor will preserve its alignment when temporarily detached during osteotomy.


New slotted plates make it easier to connect different external support (rings, threaded rods and foot plates).




Slotted Plate



TrueLok IS STABLE

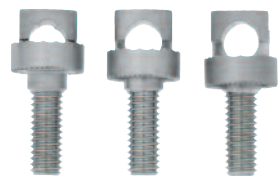

Universal Half Pin Fixation Bolt

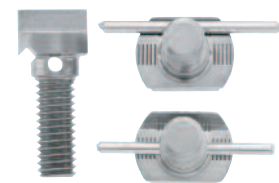

Universal Wire Fixation Bolt

The anti-rotation feature on connection elements helps prevent undesired movement.

TrueLok's unique wire and half-pin fixation elements provide efficient bone segment stabilization.

- Universal half-pin fixation bolt secures 4, 5 or 6mm half-pins.
- Universal wire fixation bolt offers cannulated and slotted options to capture wires.
- Secure connection of the wire fixation bolts dramatically reduces the need to retension wires.

A biomechanically superior metal-plastic interface is stable both during and between adjustments.



Detail of anti-rotation feature




Metal-Plastic Interface

Telescopic Linear Distractor


TrueLok IS VERSATILE

TrueLok modular external supports (rings, half-rings, 5/8 rings and arches) allow easy adjustments to the system.

- A stable, multi-planar foot support with half-rings and extension plates fit any size or shape of foot.
- The TrueLok universal lockable hinge offers self-alignment flexibility.

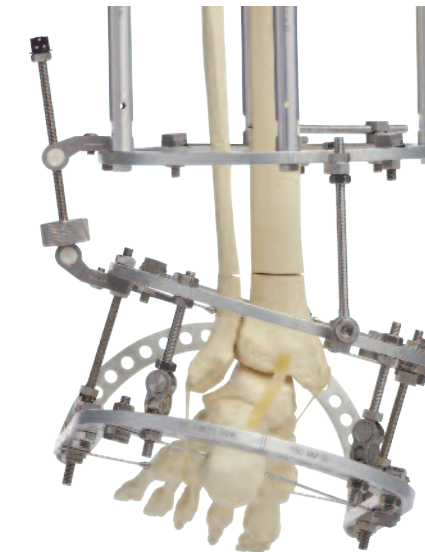

Ankle Valgus Deformity Correction

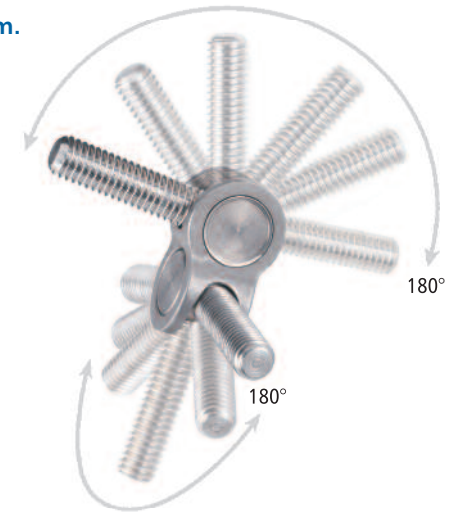

Universal Lockable Hinge

Full control of bone segments accomplished with TrueLok allows precise movement of segments without compromising stability.



Distal Femur Deformity Correction



Mid-Shaft Tibia Deformity Correction